



No. TDD001

### CellSCAFLD® 3D Cell Culture Dish

#### Description

35mm 3D Cell Culture Dish, 1 Scaffold, Surface Treated, Sterile


#### Purpose

A micro environment for cells that are similar to the in vivo conditions used in stem cells, tissue engineering, drug research and development, and cell biology

#### Materials

Dish: GPPS (General Polystyrene) 3D scaffold: GPPS (General Polystyrene)

#### Dimensions (Unit: mm)

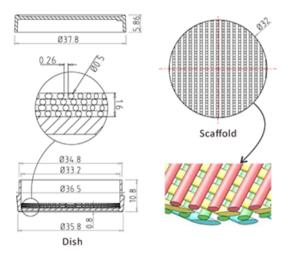

#### Features

- The scaffold is made from virgin polystyrene with a wire diameter of 500µm and a wire spacing of 260µm. It produces a large surface area than regular cell culture products and is structured with 3-dimensional channel facilitating the transmission of nutrients, consistency of metabolic activity and the accuracy of results in 3D cell culture
- Easy cell secretion collection, saving time and eliminating extra steps
- · Cytokine and growth factor resistant
- Non-autoclavable
- DNase/RNase free and Non-pyrogenic


Easy to Use

- Sterilized by irradiation SAL10<sup>-6</sup> (ISO11137)
- Shelf Life: 3 years after month of production
- Manufactured in a class 100,000 room environment
- Manufactured under ISO13485 and ISO9001 quality management system





Prepare the required volume of cell suspension

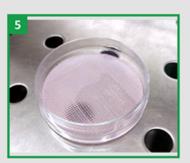


Put the dish into a 37°C and 5% CO, incubator for culturing for three hours



Add the cell suspension to the 3D scaffold slowly



After three hours, slowly add the cell culture medium through the dish internal wall

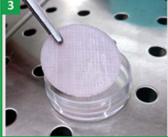

Ensure that the 3D scaffold is fully covered with cell suspension and avoid overflow



Place the 3D scaffold into the incubator once the cell culture medium covers the scaffold completely

Use tweezers to pick up the 3D scaffold and place it into the tissue culture dish

## CellSCAFLD® 3D Cell Culture Dish

Cat. No.	Туре	Fiber Diameter	Pore width	Scaffold Diameter	Scaffold Thickness	Scaffold Growth Area	Scaffold Surface Type	Sterilization	Packaging Configuration
TDD032035	1 scaffold in 35mm Dish	ø500 μm	260 μm	ø32 mm	1.6 mm	43 cm²	Treated	Yes	1/blister pack, 40/case





No. TDD002

### CellSCAFLD® 3D Cell Culture Dish

#### Description

60mm 3D Cell Culture Dish, 1 Scaffold, Surface Treated, Sterile


#### Purpose

A micro environment for cells that are similar to the in vivo conditions used in stem cells, tissue engineering, drug research and development, and cell biology

#### Materials

Dish: GPPS (General Polystyrene)
3D scaffold: GPPS (General Polystyrene)

#### Dimensions (Unit: mm)

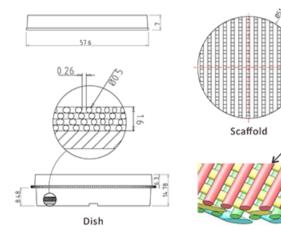

#### Features

- The scaffold is made from virgin polystyrene with a wire diameter of 500µm and a wire spacing
  of 260µm. It produces a large surface area than regular cell culture products and is structured
  with 3-dimensional channel facilitating the transmission of nutrients, consistency of metabolic
  activity and the accuracy of results in 3D cell culture
- Easy cell secretion collection, saving time and eliminating extra steps
- · Cytokine and growth factor resistant
- Non-autoclavable
- DNase/RNase free and Non-pyrogenic
- Sterilized by irradiation SAL10<sup>-6</sup> (ISO11137)
- · Shelf Life: 3 years after month of production
- Manufactured in a class 100,000 room environment
- Manufactured under ISO13485 and ISO9001 quality management system

## ander 15015 105 and 1505 001 quan

## Easy to Use




Prepare the required volume of cell suspension



Add the cell suspension to the 3D scaffold slowly



Ensure that the 3D scaffold is fully covered with cell suspension and avoid overflow

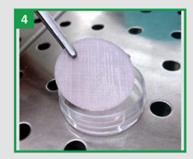

STERILE R

Use tweezers to pick up the 3D scaffold and place it into the tissue culture dish



Put the dish into a 37°C and 5% CO<sub>2</sub> incubator for culturing for three hours



After three hours, slowly add the cell culture medium through the dish internal wall



Place the 3D scaffold into the incubator once the cell culture medium covers the scaffold completely

### CellSCAFLD® 3D Cell Culture Dish

Cat. No.	Туре	Fiber Diameter	Pore width	Scaffold Diameter	Scaffold Thickness	Scaffold Growth Area	Scaffold Surface Type	Sterilization	Packaging Configuration
TDD032060	1 scaffold in 60mm Dish	ø500 μm	260 µm	ø51 mm	1.6 mm	109 cm²	Treated	Yes	1/blister pack, 30/case





No. TDD003

STERILE R

Use tweezers to pick up the 3D

scaffold and place it into the

tissue culture dish

### CellSCAFLD® 3D Cell Culture Dish

#### Description

70mm 3D Cell Culture Dish, 1 Scaffold, Surface Treated, Sterile


#### Purpose

A micro environment for cells that are similar to the in vivo conditions used in stem cells, tissue engineering, drug research and development, and cell biology

#### Materials

Dish: GPPS (General Polystyrene)
3D scaffold: GPPS (General Polystyrene)

Dimensions (Unit: mm)

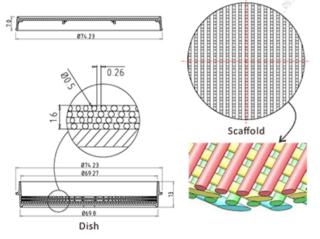


#### Features

- The scaffold is made from virgin polystyrene with a wire diameter of 500µm and a wire spacing
  of 260µm. It produces a large surface area than regular cell culture products and is structured
  with 3-dimensional channel facilitating the transmission of nutrients, consistency of metabolic
  activity and the accuracy of results in 3D cell culture
- Easy cell secretion collection, saving time and eliminating extra steps
- · Cytokine and growth factor resistant
- Non-autoclavable
- DNase/RNase free and Non-pyrogenic
- Sterilized by irradiation SAL10<sup>-6</sup> (ISO11137)
- Shelf Life: 3 years after month of production
- Manufactured in a class 100,000 room environment
- Manufactured under ISO13485 and ISO9001 quality management system

## Easy to Use



Prepare the required volume of cell suspension

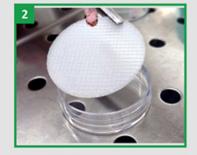


Add the cell suspension to the 3D scaffold slowly



After three hours, slowly add the cell culture medium through the dish internal wall



Ensure that the 3D scaffold is fully covered with cell suspension and avoid overflow

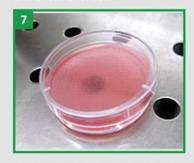

Place the 3D scaffold into the incubator once the cell culture medium covers the scaffold completely



# CO<sub>2</sub> incubator for culturing for three hours

Put the dish into a 37°C and 5%

## CellSCAFLD® 3D Cell Culture Dish

Cat. No.	Туре	Fiber Diameter	Pore width	Scaffold Diameter	Scaffold Thickness	Scaffold Growth Area	Scaffold Surface Type	Sterilization	Packaging Configuration
TDD032070	1 scaffold in 70mm Dish	ø500 μm	260 µm	ø67.5 mm	1.6 mm	191 cm²	Treated	Yes	1/blister pack, 30/case